



Second Opinion Professional

How to Manually send a Second Opinion Folder

To manually send a Second Opinion Folder (.SOD File), select the content to be sent, then click the Export Folder choice located in the Folder menu.

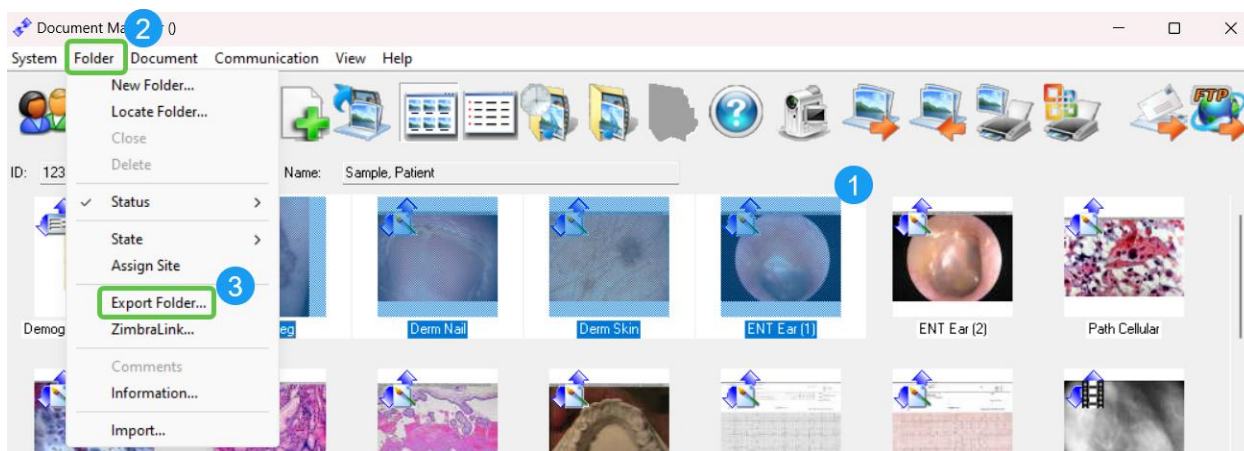


The Export dialog box that comes up is the same one that is launched when using FTPSend.

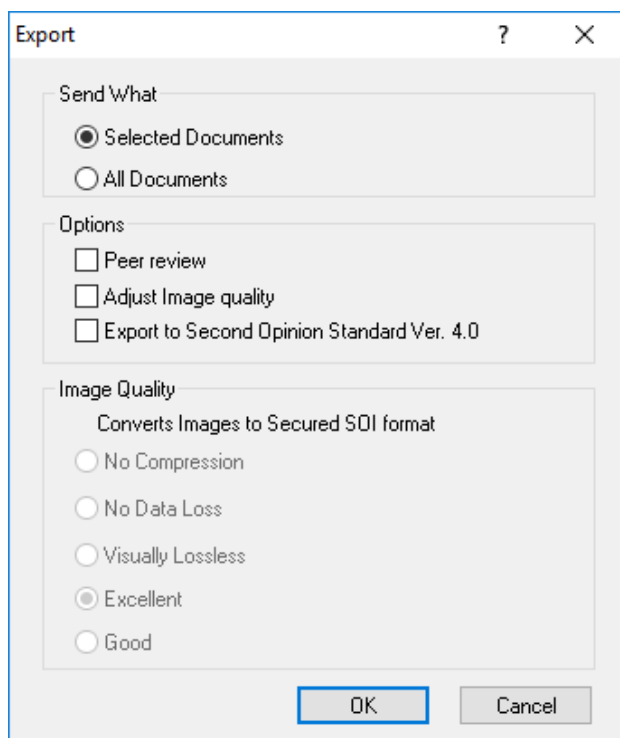



When finished, name the SOD file and save it in the desired location.

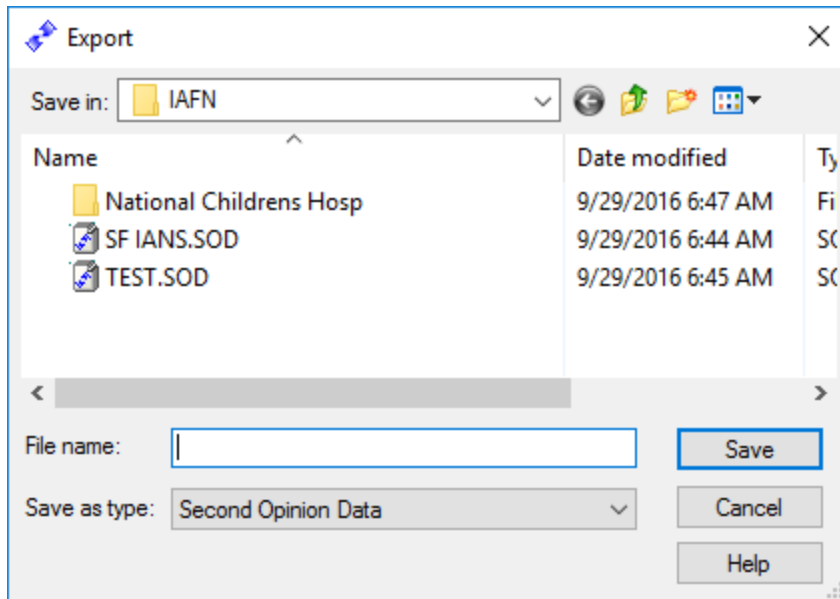

This encrypted .SOD file can now be attached to an email or saved to a USB drive. The recipient of the .SOD file must have the Second Opinion Professional software to open the file.

For More Information or Technical Support, Give Us a Call

Second Opinion Telemedicine Solutions, Inc.

310-802-6300

www.sotelemed.com

corp@sotelemed.com

