



Second Opinion Professional

How to generate a PDF document for auto import into a Patient's Chart

This document explains the steps necessary to configure the WReport feature so it will export a PDF document and save it to a configured location. From there it can be auto imported into a Patient's medical record.

Configure Word Report to Export

Word Report (WReport) is a program that programmatically interacts with Microsoft Word to merge Second Opinion images and data into specialized Word templates that were created using the Second Opinion Report Template Designer.

To activate the feature, using notepad or some other text editor, create a WReportCharting.ini file and insert the following text.

[Charting]

Location=<insert auto import folder location>

Below are a couple of examples to use for the Location field.

Note: UNC Server paths are preferred but it also supports Mapped Drive locations.

Location=\\Someserver\folder\import

OR

Location=C:\Users\Public\Second Opinion\temp

When ready, place the completed file into one of two locations.



C:\Users\Public\Second Opinion\WreportCharting.ini (Configures a single Workstation)

Or


The Second Opinion Document folder. This location can be found by looking at the DocumentPath= contained in the setting.ini file.

C:\Users\Public\Second Opinion\Document\WreportCharting.ini (Configures all network Workstations)

To test it, create a document using a Word Report Template. When ready, exit Word. Insert the Report into the Patient's folder and lastly click Yes to Save the document into the Patient's Chart.

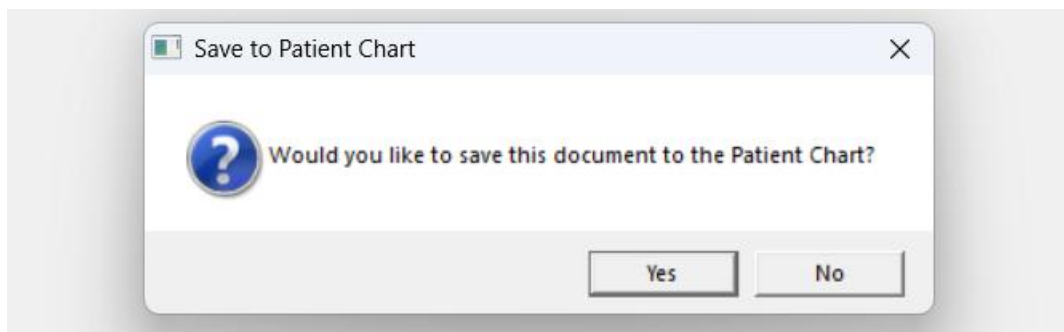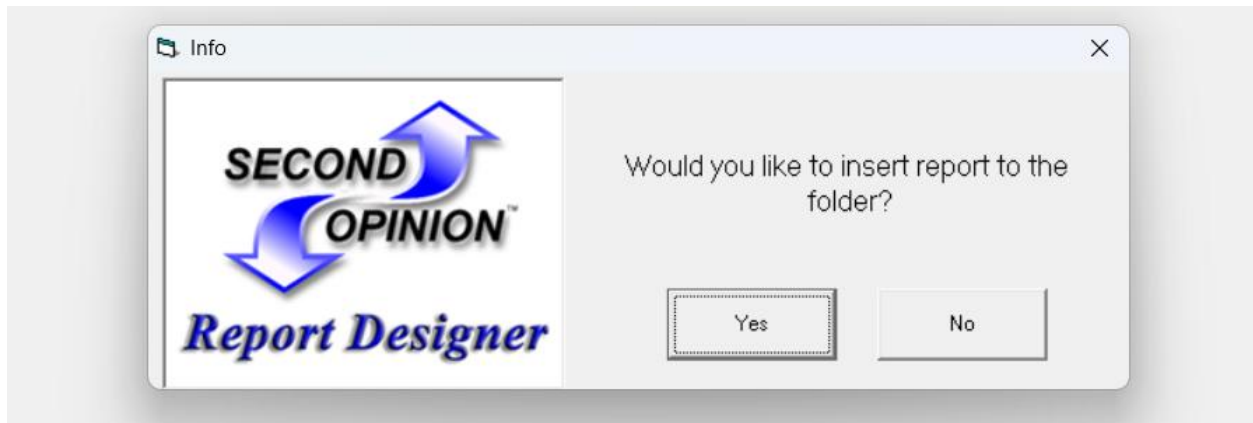

When finished, go to the configured location and verify the PDF report file is there.

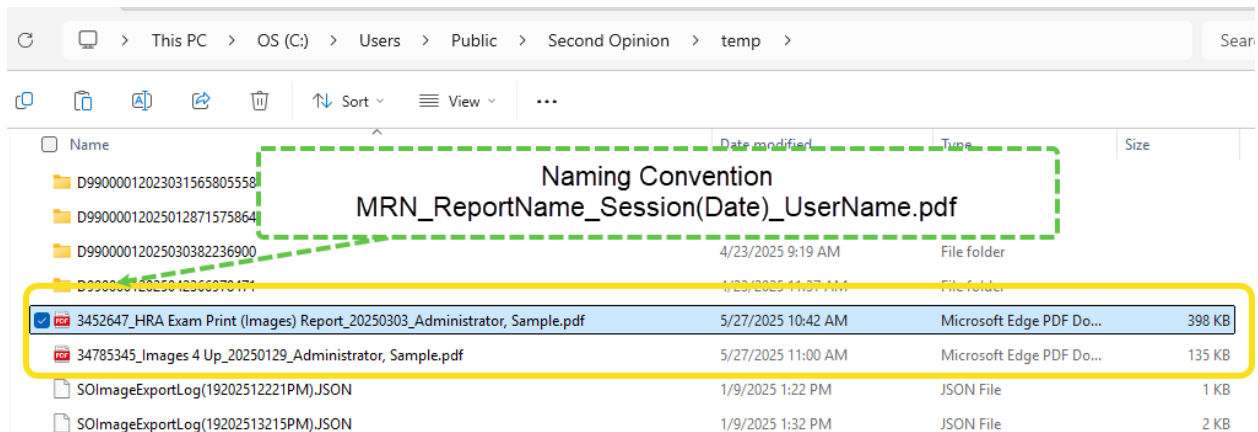



(This Page is Intentionally Left Blank)

For More Information or Technical Support, Give Us a Call

Second Opinion Telemedicine Solutions, Inc.

310-802-6300

www.sotelemed.com

corp@sotelemed.com

